

Fast + Simple
Focused on Veterinary Diagnostics

FASTest® FeLV-FIV ad us. vet.

World-wide most common virus infection in cats

Fast test for the qualitative detection of Feline Leukaemia Virus (FeLV) antigens and/or Feline Immunodeficiency Virus (FIV) antibodies in whole blood, plasma or serum of the cat

Fast and reliable diagnostics

FeLV: aetiological antigen detection e.g. immunosuppression, all kinds of infections, tumours

FIV: antibody detection e.g. immunosuppression, lymphadenitis, all kinds of chronic or recurrent inflammations, anaemia

FeLV/FIV exclusion diagnostics
in case of unspecific clinical symptoms

Annual check-up

- Routine test before FeLV/FIV vaccination
- Determination of FeLV/FIV status

Routine test of risk groups

- New additions of unknown origin
- Uncontrolled multiple cat households
- Stray cats



- Simple test procedure with whole blood, plasma or serum
- Fast test interpretation after 10 minutes
- Reliable clinical diagnostics
- FeLV: Sensitivity 95.0% & Specificity 99.9%
- FIV: Sensitivity 96.4% & Specificity 99.2%
- Storage at room temperature (15–25 °C)
- Long shelf life
- Compact test box with 2, 6, 15, 25 or 50 tests

FASTest[®] FeLV-FIV ad us. vet.

The Feline Leukaemia Virus (FeLV) as well as the Feline Immunodeficiency Virus (FIV) occur world-wide in felids. Prevalences in Germany range from 2% (FeLV) to 2–3% (FIV).

Both diseases cause similar clinical symptoms in infected cats due to their immunosuppressive effects. FeLV antigen or FIV antibody detection are proven to be the diagnostic method of choice confirming a FeLV or FIV infection, respectively.


FeLV infected, viraemic cats normally show high concentrations of extracellular (free) FeLV antigen 3 weeks post infection. About 95% of the FIV infected cats show high FIV antibody concentrations approximately 4 weeks post infection.

Based on highly specific monoclonal FeLV antibodies and recombinant FIV proteins, respectively, **FASTest[®] FeLV-FIV** is an important diagnostic tool for the evaluation of clinically suspicious cats.

Test procedure



Test interpretation

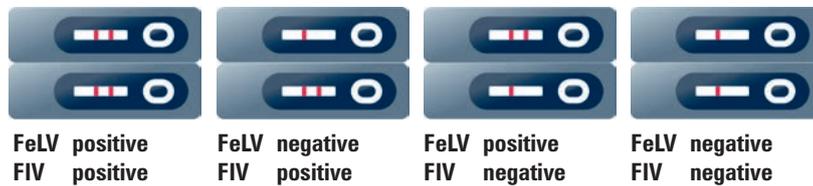

FASTest[®] FeLV is recommended in case of single testing of FeLV before vaccination.

FASTest[®] FIV is recommended in case of single testing of FIV before vaccination.

Distribution:

EN 10-2019

